



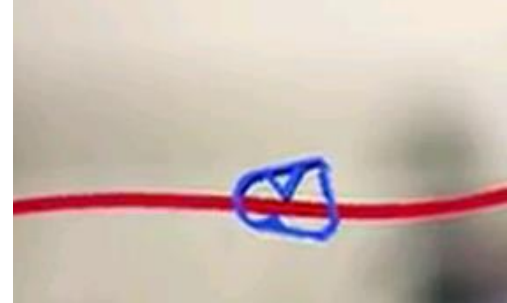




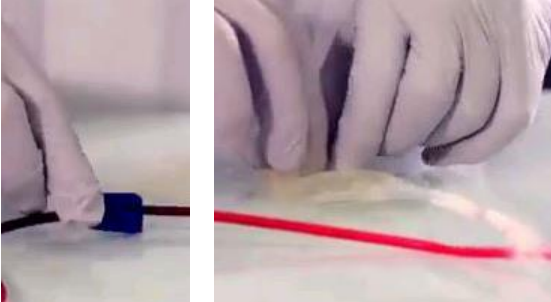
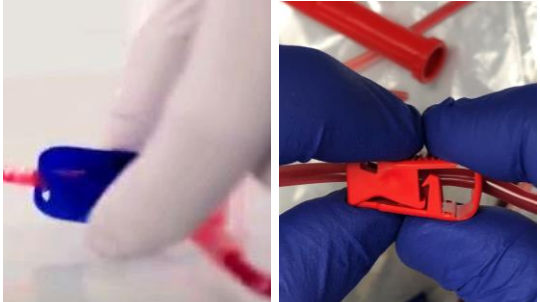



Pour copie papier seulement : Le document original est approuvé avec signature numérique sous la responsabilité de CFA.

<p>1 Ensure the mother does not present any signs of infection or evidence of illegal/unprescribed drug use</p> 	<p>2 Close the 2 blue clamps</p> 	<p>3 Disinfect the cord unidirectionally</p> 
<p>4 Insert the needle into the vein of the cord</p> 	<p>5 Open the 2 blue clamps and lower the collection bag</p> 	<p>6 After the blood has been collected, milk the cord</p> 
<p>7 Close the 2 blue clamps prior to removing the needle</p> 	<p>8 Push up and lock the needle protector</p> 	<p>9 Break the cannula from the anticoagulant pouch</p> 
<p>10 Open the proximal blue clamp and roll the anticoagulant pouch up onto itself</p> 	<p>11 Close the proximal blue clamp and the permanent red clamp</p> 	<p>12 Sign the collection form</p> 

Reference: PFN-00028

LIST OF MODIFICATIONS:

SECTION	MODIFICATIONS	JUSTIFICATION	ORIGIN OF INFORMATION
11	Replacement of the red clamp picture	Update following its replacement	CC-PRB-01341
Ref. CFA : OC-06161			